

[root](#) > [In practice](#) > [How it's done](#) >




Exact saw-cut models produced quickly and simply.

**8 minutes working time
means economy and efficiency:**



The impression aligned in the model-positioner and fixed with adhesive plasticine is cast in a single process using the **model-tray** base former.

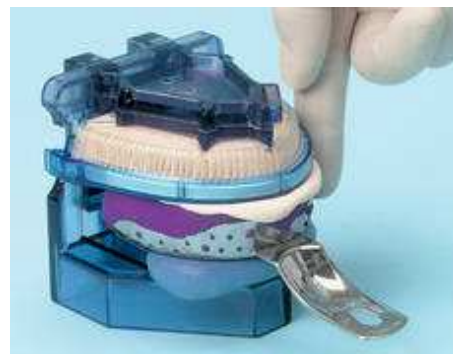

After 30 minutes*, lift the impression off the model ...

... after 90 minutes*, remove the model from the mould



... fast, precise sawing with the **model-cut** sawing machines ...



- wafer-thin saw cuts
- precise and strong cut segments thanks to asymmetric guide grooves
- the palatal part is preserved



* from the start of mixing the plaster

Julius-Vosseler-Str. 42 | DE-22527 Hamburg | Telefon 040-3990366-0 | Fax 040-3990366-22 | info@model-tray.de | [Impressum](#)